



## Hydraulic Medical Chair

OC-1A  
OC-1B

Aiming at new levels in quality





## Ideal chair for operators

### Hydraulic Medical Chair

# OC-1A/OC-1B

#### Hard and thick cushion

Since hard and thick cushion is used, the body does not sink very deep and the operator doesn't feel so tired even after a long time sitting.

#### Special seat shape

The seat is of special triangle which reduces the pressure on the rear thighs.

#### Backrest to support the waist


The backrest, with a width and curvature designed to hold the upper part of the waist, not only supports the back but also enables the operator to easily maintain a posture that lightens the load on the waist.

#### Simple elevating section

Since the existing simple, highly durable mechanism is adopted to the pedal type elevating section, raising, lowering and seat rotation lock can be performed easily and securely.

#### Armrest with cushion

The OC-1A armrest with cushion is slightly inclined inwards so that the arms are placed naturally.

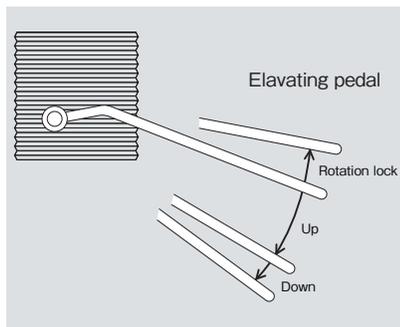


#### Caster with big diameter and wide width

Casters with big diameter and wide width reduce the resistance and enable the operator to move easily.

#### Stable base section

The base is designed to be stable with five arm legs and the hydraulic mechanism of low center of gravity.

### Dimensions

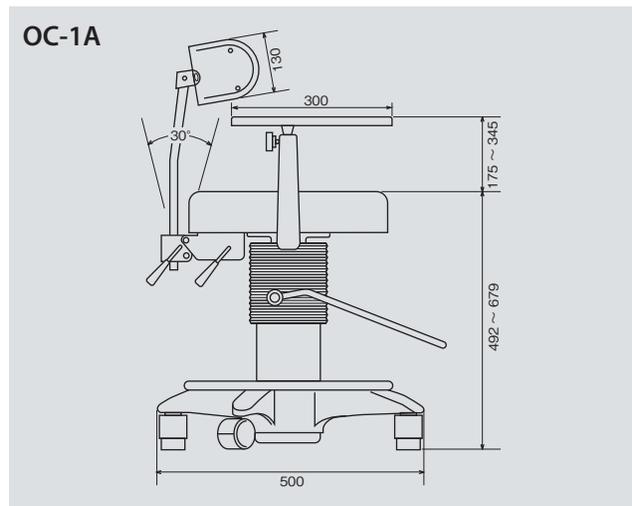

### Specifications

Elevation method	Pedal, hydraulic type
Elevation stroke	187mm
Backrest up/down stroke	230mm
Seat size	380mm(W)×390mm(D)
Armrest front/rear stroke	OC-1A : 170° (250mm) OC-1B : Not provided
Weight	OC-1A : 31kg OC-1B : 25.5kg

● Design and specifications are subject to change as improvements are made to the product.

**TAKAGI SEIKO CO.,LTD.**

330-2 IWAFUNE, NAKANO-SHI, NAGANO-KEN, 383-8585, JAPAN  
TEL.+81-269-22-4512 FAX.+81-269-26-6321  
URL:<http://www.takagi-j.com> E-mail:[info@takagi-j.com](mailto:info@takagi-j.com)

